# Masimo SafetyNet Alert™

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## For Export Only, Not for Sale in the USA

There may be information provided in this manual that is not relevant for your system. Do not operate Masimo SafetyNet Alert without completely reading and understanding these instructions. If you encounter any serious incident with product, please notify the competent authority in your country and the manufacturer.

Wireless Radio:

Contains FCC ID: 2AC7Z-ESP32SOLO1; FCC ID: VKF-CONNHUB

Contains IC: 21098-ESP32SOLO1; IC: 7362A-CONNHUB

Masimo Corporation  
52 Discovery  
Irvine, CA 92618, USA  
Tel.: 949-297-7000  
Fax.: 949-297-7001  
www.masimo.com

EU authorized representative for Masimo Corporation:

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Patents: [www.masimo.com/patents.htm](http://www.masimo.com/patents.htm)

Masimo® and  are federally registered trademarks of Masimo Corporation.

Masimo SafetyNet™ and Masimo SafetyNet Alert™ are trademarks of Masimo Corporation.

All other trademarks and registered trademarks are property of their respective owners.

© 2022 Masimo Corporation

# Contents

---

About This Manual.....	5
Getting Started Guides.....	5
In-App Videos.....	5
Product Description and Intended Use.....	5
Product Description.....	5
Who the Device is intended to be used on (Indications for Use).....	6
Intended Use.....	6
Indications for Use.....	6
Contraindications.....	6
Safety Information.....	7
Safety Warnings.....	7
Masimo Sensor.....	8
Performance Warnings.....	8
Masimo Sensor.....	10
Cleaning and Service Warnings.....	10
Compliance Warnings.....	10
Risks and Benefits.....	12
Your Body and Oxygen.....	12
Risks of Masimo SafetyNet Alert.....	13
Benefits of Masimo SafetyNet Alert.....	14
Description.....	15
Masimo SafetyNet Alert System.....	15
Features.....	16

Home Medical Hub Overview .....	16
Masimo Sensor and Masimo Chip Overview.....	18
Basic Setup and Use.....	19
Getting Started.....	19
Step 1: Smart Phone Setup .....	19
Step 2: Masimo SafetyNet Alert System Setup.....	21
Home Medical Hub Setup .....	21
Masimo Sensor Setup.....	22
Add Emergency Contacts .....	28
Step 3: View Live Data.....	29
Alarms and Notifications.....	30
Advanced Use .....	35
Health Settings .....	35
Change Option .....	35
Oxygen Level.....	36
Troubleshooting.....	39
Masimo SafetyNet Alert Messages.....	39
Troubleshooting Masimo SafetyNet Alert.....	45
Troubleshooting Measurements.....	50
Appendix.....	51
Specifications.....	51
Masimo SafetyNet Alert App .....	51
Masimo Sensor.....	51
Masimo Chip .....	54
Home Medical Hub .....	54

Symbols.....	63
Citations.....	64
Service and Maintenance.....	65
Cleaning.....	65
Customer Support.....	66
Concepts of Alarm Response Delay.....	67
Index/Indice.....	143

## About This Manual

---

### Getting Started Guides

Use the quick reference guide included for the following:

- Basic Setup and starting your first monitoring session.

### In-App Videos

Watch the videos in your app to find out more about:

- Setting up your Masimo SafetyNet Alert System.
- Masimo Sensor placement.
- Connecting your Masimo SafetyNet Alert System to WiFi.

## Product Description and Intended Use

---

### Product Description

**Masimo SafetyNet Alert** - a system that provides spot-checking and continuous monitoring of physiological data by wireless communication to medical technologies. The system includes the following pieces:

**Masimo SafetyNet Alert App** - software application installed on a smart phone that provides the graphical user interface to display your data and alarm condition status.

**Masimo Sensor** - wireless, wearable sensor that provides the physiological data.

**Masimo Chip** - attaches to the wireless Masimo Sensor and connects wirelessly to the Home Medical Hub.

**Home Medical Hub** - device that communicates monitoring data wirelessly from the medical technologies to the Masimo Cloud.

**Masimo Cloud** - a server accessed over the internet that gathers and stores measured data communicated wirelessly from the Home Medical Hub.

## Who the Device is intended to be used on (Indications for Use)

### Intended Use

The Masimo SafetyNet Alert is intended for the continuous monitoring of physiological parameters through communication to compatible medical technologies. Collected data from compatible sensors are transferred and displayed via the subject device, providing a system that allows physicians to monitor levels of SpO<sub>2</sub> desaturation notifications.

### Indications for Use

The Masimo SafetyNet Alert is intended for the continuous monitoring of physiological parameters through communication to compatible medical technologies.

The Masimo SafetyNet Alert is indicated for the continuous monitoring of functional oxygen saturation of arterial hemoglobin (SpO<sub>2</sub>) and Pulse Rate (PR) in adult, pediatric, infant, and neonates during motion, no-motion and low perfusion conditions in hospital and home environments.

The Masimo SafetyNet Alert is also indicated for the monitoring of adults and pediatrics for substance induced oxygen desaturations (e.g. OIRD caused by oral or injectable opioids) in hospital and home environments.

### Contraindications

Masimo SafetyNet Alert is contraindicated for the monitoring of opioid users on supplemental oxygen greater than 2.0 L/min.

# Safety Information

---

Before using the Masimo SafetyNet Alert, read the following safety information carefully.

## Safety Warnings

**WARNING:** For safe use, do not use if it appears damaged.

**WARNING:** For your safety, do not try to repair a device that is damaged.

**WARNING:** Place the Home Medical Hub where it will not fall on anyone.

**WARNING:** Make sure the oxygen level is displayed to ensure that monitoring has started.

**WARNING:** For safe use avoid the following:

- Do not place the Home Medical Hub on a wet surface.
- Do not soak any part of the system in liquid.
- Do not try to sterilize.
- Only use Masimo recommended solutions for cleaning your system. See **Cleaning** on page 65.
- Do not clean while the device is in use.

**WARNING:** Keep small parts away from small children. Small items can be a choking hazard.

**WARNING:** Carefully position any cables to avoid possible strangulation or entanglement.

**WARNING:** Do not place the Home Medical Hub in areas where there are a lot of flammable gases such as anesthetics, oxygen, or nitrous oxide present to prevent risk of fire.

**WARNING:** During defibrillation, avoid contact with any devices applied to the person to prevent electric shock.

**CAUTION:** Place the Home Medical Hub where you can easily disconnect it from AC power in case of an emergency.

**CAUTION:** Only use the AC power supply included with your Home Medical Hub to prevent damage to the device.

**Note:** Do not monitor more than one person at a time with Masimo SafetyNet Alert.

# Masimo Sensor

**WARNING:** Avoid wrapping the sensor too tightly around your finger to avoid injury.

**WARNING:** To avoid skin injury, consider moving the sensor to another finger after it has used more than 8 hrs.

## Performance Warnings

**WARNING:** Do not self-diagnose or self-medicate on the basis of the measurements. Always consult your physician.

**WARNING:** Do not use Masimo SafetyNet Alert as an apnea monitor. Masimo SafetyNet Alert is not designed to detect apnea.

**WARNING:** Only use Masimo-approved parts with Masimo SafetyNet Alert to make sure the device works correctly.

**WARNING:** Check your device is working before use for longer periods by checking to see your starting oxygen level.

**WARNING:** Avoid placing liquids on or near the Home Medical Hub. Liquids spilled on Masimo SafetyNet Alert may affect how it works.

**WARNING:** Do not place anything on top of the Home Medical Hub to avoid damaging or blocking the wireless signal or muffling the audible alarms.

**WARNING:** The following factors can cause false oxygen measurements:

- Sensor is applied incorrectly.
- Restricted blood flow to the sensor site.
- Sensor site has nail polish, acrylic nails, glitter, etc.
- Sensor site has moisture, birthmarks, skin discoloration, or foreign objects.
- Placed near other devices that may interfere with its operation.
- Excessive movement of the sensor is present.
- Possible medical conditions that might interfere with readings.

**WARNING:** To ensure proper notification function, check the following occasionally:

- Notification features is turned on your Smart device (i.e., sounds, vibrations, etc.).

- Smart device battery is fully charged or plugged in.
- Oxygen values are displayed on the App live view.

**CAUTION:** Do not place your smart phone where someone can change the App settings without you knowing.

**CAUTION:** Avoid using Masimo SafetyNet Alert under bright light sources and direct sunlight to maintain the performance of the device.

**CAUTION:** Check that the Masimo SafetyNet Alert alarm can be heard from other rooms in your home, especially when noisy appliances such as vacuum cleaners, dishwashers, clothes dryers, televisions, or radios are operating.

**CAUTION:** Keep the Home Medical Hub plugged in while in use. Loss of power may limit the notifications available.

**CAUTION:** Only use the AC power supply that came with the Home Medical Hub to prevent damage to the device.

**CAUTION:** Do not connect to an electrical outlet controlled by a wall switch or dimmer.

**CAUTION:** When using Masimo SafetyNet Alert, locate the devices away from sources that may interfere with the Bluetooth connection. The presence of other devices that may create radio frequency interference (RFI) may result in loss of Quality of Service of the Bluetooth connection (see Specifications for details). Devices that may cause RFI include but are not limited to the following: cell phones, laptops and tablets, pagers, Bluetooth devices, devices with remote controls, electrocautery equipment, diathermy equipment, and baby monitors.

**CAUTION:** To ensure security and prevent tampering of your smart phone, while using the Masimo SafetyNet Alert:

- Smart device should be located with responsible users.
- Smart device should not be left unattended.
- Security features on smart device should be activated.
- Unauthorized changes should not be made to the Masimo SafetyNet Alert system.

**CAUTION:** Check your setup by viewing the display on the App and the Home Medical Hub. The App and Home Medical Hub will provide an indication if there is a problem with the internet connection.

**CAUTION:** Check your system setup by viewing the display on the smartphone Masimo SafetyNet Alert Application. The Masimo SafetyNet Alert App will provide an indication if there is a problem with the internet connection.

**CAUTION:** Avoid placing Masimo SafetyNet Alert near other wireless devices to prevent loss of performance.

## Masimo Sensor

**WARNING:** Properly apply the Masimo Sensor according to sensor's Directions for Use. Applying the sensor incorrectly could result in incorrect or no readings.

**CAUTION:** Avoid using Masimo Sensor under bright light sources and direct sunlight to maintain the performance of the device.

**CAUTION:** To maintain Bluetooth connectivity with the Masimo Sensor, ensure that the Home Medical Hub is within the specified distance and in line-of-sight of the Masimo Sensor. See *Specifications* on page 51.

## Cleaning and Service Warnings

**WARNING:** Do not attempt to remanufacture, recondition or recycle the Masimo Sensor, Masimo Chip or Home Medical Hub to prevent harm or damage to the system.

**WARNING:** Always turn off and unplug the Home Medical Hub before cleaning to prevent harm or damage to the device.

**CAUTION:** Do not clean Masimo Sensor, Masimo Chip or Home Medical Hub with undiluted bleach, petroleum-based products, acetone, or other harsh solvents. Clean only with the solutions specified in this manual to prevent damage to the device.

**CAUTION:** Do not submerge Masimo Sensor, Masimo Chip or Home Medical Hub in liquid or attempt to sterilize by any method to prevent damage to the device.

## Compliance Warnings

**WARNING:** Changes or modifications not approved by Masimo can void the user's authority to operate the equipment.

**WARNING:** The frequency bands of this device (2.4 GHz and 5.15 to 5.25 GHz) are only for indoor use in accordance with international telecommunication requirements.

**WARNING:** Only use Masimo authorized devices with Masimo SafetyNet Alert. Using unauthorized devices with Masimo SafetyNet Alert may result in damage to the device and/or patient injury.

**WARNING:** Portable RF communications equipment (including peripherals such as antenna cables and external antennas) should be used no closer than 30 cm (12 inches) to any part of the Masimo SafetyNet Alert, including cables specified by the manufacturer. Otherwise, degradation of the performance of this equipment could result.

**CAUTION:** Disposal of product: Comply with local laws when disposing the device and/or its accessories.

**CAUTION:** Do not place the Masimo SafetyNet Alert near electrical equipment that may affect the device, preventing it from working properly.

**CAUTION:** To minimize radio interference, other electrical equipment that emits radio frequency transmissions should not be near Masimo SafetyNet Alert.

**CAUTION:** Keep the Masimo SafetyNet Alert away from electrical equipment that emits radio frequencies to minimize radio interference. Radio interference may result in no or inaccurate readings.

**Note:** Masimo SafetyNet Alert complies with the limits for a Class B digital device, per Part 15 of the FCC Rules. These limits were designed to provide reasonable protection against harmful interference in a residential installation. Masimo SafetyNet Alert generates, uses, and can radiate radio frequency energy and may cause interference with radio communications. To determine if Masimo SafetyNet Alert interferes with radio or television reception, turn it off and see if the interference stops. To correct the interference, try the following:

- Adjust or move the receiver's antenna.
- Move the receiver farther away from Masimo SafetyNet Alert.
- Plug the receiver and Masimo SafetyNet Alert into outlets on different circuits.
- Consult the dealer or a radio/TV technician for help.

**Note:** This equipment has been tested and found to comply with the Class B limits for medical devices according to the IEC 60601-1-2: 2014. These limits are designed to provide reasonable protection against harmful interference in all establishments, including domestic establishments.

**Note:** This device complies with part 15 of the FCC Rules and Industry Canada’s license-exempt RSS standards. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

**Note:** To satisfy RF exposure requirements, this device and its antenna must operate with a separation distance of at least 20 cm from all persons and must not be co-located or operating in conjunction with any other antenna or transmitter.

**Note:** This device complies with Industry Canada license-exempt RSS standard(s). Operation is subject to the following two conditions: (1) this device may not cause interference, and (2) this device must accept any interference, including interference that may cause undesired operation of the device.

**Note:** Users are advised that high-power radars are allocated as primary users (i.e., priority users) of the bands 5250-5350 MHz and 5650-5850 MHz and that these radars could cause interference and/or damage to LE-LAN devices.

**Note:** When using Masimo SafetyNet Alert consideration should be taken to local government frequency allocations and technical parameters to minimize the possibility of interference to/from other wireless devices.

## Risks and Benefits

---

### Your Body and Oxygen

When we breathe, oxygen moves in our blood around our body. When we are healthy or not fatigued, our lungs are able to bring the oxygen needed by our brain, organs, and tissue. As the state of our body is always changing, it is common to see continual fluctuations in your blood oxygen over time.

However, when we are sick, fatigued, or having trouble breathing, the body can require more oxygen than can be delivered. When this occurs, your blood oxygen level can drop. These changes cannot be seen, but if our oxygen gets too low, you may feel the effects such as being lightheaded, dizzy, or having difficulty breathing. If your oxygen level stays too low for too long, the lack of oxygen to your brain, organs, and tissue or hypoxia can have serious long-term effects to your health including death.

The Masimo SafetyNet Alert allows you to see your oxygen level so that you can understand how your body and oxygen are connected.

If you have a medical condition that affects your breathing (e.g., asthma, COPD), your normal blood oxygen level may be lower. For these conditions, it is important that you seek the advice of your physician.

For other conditions where your blood oxygen becomes low, it is good to remember that our blood oxygen is connected to our breathing and always changing.

If you have any concerns that your oxygen level is constantly decreasing, staying low, or you are not feeling well, it is very important to not self-diagnose and seek the advice of your physician.

When using any medical device, there are risks and benefits.

## Risks of Masimo SafetyNet Alert

- As with all medical electrical devices, there is always a risk of electrical, mechanical, and fire hazards. However, these risks have been mitigated through the design and testing of Masimo SafetyNet Alert.
- As with all devices with small parts, there is always a risk of a child's swallowing a component or choking. Keep small parts away from small children to prevent swallowing or choking.
- Items applied directly to the skin may cause skin irritation, pressure injury, or general discomfort. Periodically check the area where the Masimo sensor is applied to prevent potential irritations.
- As with all types of alarms, there is always a risk of missing or not hearing an alarm. To minimize this risk, ensure Masimo SafetyNet Alert is placed where it can be heard and the speaker is not blocked by other objects.
- As with all types of alarms, there is always a risk of false alarms. To minimize false alarms, ensure you follow all device and Masimo sensor directions, and follow troubleshooting instructions. Never assume an alarm is false; promptly address every alarm.

## Benefits of Masimo SafetyNet Alert

- Masimo SafetyNet Alert has the ability to monitor wirelessly without being tethered to a monitor. This allows users the freedom of movement during activities, including sleep, not conveniently possible with a cabled connection.
- Masimo SafetyNet Alert is equipped with audible and visual physiological alarms that can alert you to changes in the user's condition while monitoring within a home environment.
- Masimo SafetyNet Alert uses the same proven pulse oximetry used in hospitals shown to have fewer false alarms<sup>1</sup> compared to other brand technology.
- Masimo SafetyNet Alert allows for the storage and review of your blood oxygen level as you go about your normal activities, including exercise or sleep.

<sup>1</sup> *Barker S.J. Anesth Analg. 2002 Oct;95(4):967-72.*

# Description

---

## Masimo SafetyNet Alert System

The Masimo SafetyNet Alert system consists of the following:

Item	Description
1	Home Medical Hub
2	Smart phone with Masimo SafetyNet Alert App Installed *
3	Masimo Sensor
4	Masimo Chip

\* Smart phone not included.

# Features

The following describes the Masimo SafetyNet Alert system parts.

## Home Medical Hub Overview

The Home Medical Hub is a device that transfers measured data from the wireless Masimo Sensor to the Masimo Cloud. It transfers data from the wireless Masimo Sensor to the Masimo Cloud and provides audible and visible alarms.

### 1. Pairing Symbol

Location on the Home Medical Hub for pairing the Masimo Chip.

### 2. Home Medical Hub Button

Used for pairing and silencing alarms.

### 3. Bluetooth Status LED

Shows the Home Medical Hub Bluetooth status. See *Home Medical Hub Lights* on page 17.

### 4. Home Medical Hub Status LED

Shows the Home Medical Hub power and monitoring status. See *Home Medical Hub Lights* on page 17.

### 5. Wi-Fi Status LED

Shows the Home Medical Hub Wi-Fi status. See *Home Medical Hub Lights* on page 17.

### 6. USB Power Connector

Power is provided to the Home Medical Hub from the AC adapter and a USB cable.

### 7. Masimo Chip Storage Tray

Location on the Home Medical Hub for storage of the chip when not in use.

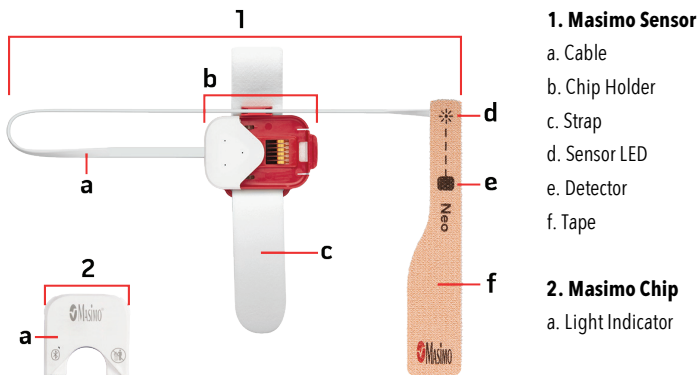
# Home Medical Hub Lights

The LED lights on the Home Medical Hub show the status of Bluetooth pairing and Wi-Fi connections.

Indicator/LED Color	What does it mean?
<b>Home Medical Hub Status LED</b>	
Solid White	Home Medical Hub is on and needs to be set up.
Solid Green	Home Medical Hub is paired with the sensor, connected to the Wi-Fi and communicating with the Masimo Cloud.
Flashing Orange	Active Level 1 Alert. See <i>Alarms and Notifications</i> on page 30.
Solid Orange	Level 1 Alert is acknowledged.
Flashing Red	Active Level 2 or 3 Alert. See <i>Alarms and Notifications</i> on page 30.
Solid Red	Level 2 or 3 Alert is acknowledged.
Flashing Orange	Sensor battery ports are blocked.
Flashing Red	A Home Medical Hub fault has been detected.
<b>Bluetooth Status LED</b>	
Solid White	Home Medical Hub is on and needs to be set up.
Off	Home Medical Hub is paired with the sensor.
Flashing White	Home Medical Hub is pairing with the sensor.
Solid Orange	Sensor is disconnected from the Home Medical Hub.
<b>Wi-Fi Status LED</b>	
Solid White	Home Medical Hub is on and needs to be set up.
Off	Home Medical Hub is connected to the Wi-Fi network.
Flashing White	Home Medical Hub is searching for or connecting to the Wi-Fi network.
Solid Orange	Home Medical Hub is connected to the Wi-Fi and the Sensor is no longer connected.

# Masimo Sensor and Masimo Chip Overview

The Masimo Sensor and Masimo Chip are for use with the Masimo SafetyNet Alert system.

## Sensor Lights

When inserted in the sensor, the light indicator on the Masimo Chip shows the status of the sensor.

Light Indicator Color	What does it mean?
Flashing Blue	Sensor is paired with the Home Medical Hub.
Flashing Purple	Sensor battery ports are blocked.
Flashing Green	Sensor is not paired with the Home Medical Hub.
Flashing Orange	Sensor battery is low.
Flashing Red	Sensor battery is very low.

# Basic Setup and Use

---

## Getting Started

The Masimo SafetyNet Alert system can be ready to use in the following steps:

1. Setup the Smart Phone - Prepare the smart phone for use with the Masimo SafetyNet Alert App.
  - Download the Masimo SafetyNet Alert App on your smart phone and complete the registration and login process.
  - Select the option that applies to the user for the Masimo SafetyNet Alert App. The option determines when alerts and notifications are triggered.
2. Setup the Masimo SafetyNet Alert System – Use the Masimo SafetyNet Alert App and complete the following:
  - Connect the Home Medical Hub using local Wi-Fi with the Masimo Cloud.
  - Pair the Masimo Sensor with the Home Medical Hub.
  - Select emergency contacts.
3. View Live Data – The Masimo SafetyNet Alert App and phone are ready for monitoring.

## Step 1: Smart Phone Setup

### Prepare the Smart Phone for Use

A compatible smart phone is required to install and operate the Masimo SafetyNet Alert App.

To use your smart phone with the Masimo SafetyNet Alert App check the following:

- **Compatibility**

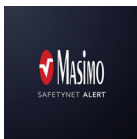
**Note:** For a list of smart phones and operating systems that work with the Masimo SafetyNet Alert App, check [www.masimo.ca/support](http://www.masimo.ca/support) before upgrading the smart phone or its operating system.

- **Smart Phone Battery is charged**
- **Bluetooth is ON**

- **Wi-Fi is ON and the smart phone is connected to the internet**
- **Time is set to the Current Local Time**

### **Download and install the Masimo SafetyNet Alert App on the smart phone:**

Download and install the Masimo SafetyNet Alert App. Follow the on-screen instructions to install. For more on how to install an app, see the smart phone's manual.

**Note:** If the Masimo SafetyNet Alert App requests for the smart phone to share its location, select *OK* or *Allow*. The smart phone location is required for Bluetooth connection.

1. After installation, open the Masimo SafetyNet Alert App and login to the user account
  - If you do not have an existing account, select **Sign Up**. Follow the on-screen instructions to create a new account and login.
2. The app walks you through the steps to set up Masimo SafetyNet Alert starting with the Home Medical Hub.

**Note:** Make sure your smartphone Bluetooth connection is turned on to connect to Masimo Sensor. Refer to the smart phone's instructions to change its settings.

### **Select an Option:**

After installing the App and setting up your account, select an Option that best applies to you. Available alert options will change based upon your selection. For complete information see ***Change Option*** on page 35. Available options are:

**COVID-19** - Select if you are recovering from COVID-19.

**Opioids** - Select if you are taking opioids.

**Other** - Select if you are not recovering from COVID-19 or not taking Opioids.

## Step 2: Masimo SafetyNet Alert System Setup

After the Masimo SafetyNet Alert App is installed and setup, follow the app instructions to set up your Masimo SafetyNet Alert System:

- Pair the App to the Home Medical Hub.
- Connect the Home Medical Hub to the Wi-Fi.
- Setup the Masimo Sensor and pair to the Home Medical Hub.
- Add your Emergency Contacts.

### Home Medical Hub Setup

Follow the instructions below to set up the Masimo SafetyNet Alert system. The Home Medical Hub is set up in two (2) steps.

**Step 1:** Pair the Masimo SafetyNet Alert App to the Home Medical Hub using Bluetooth.

**Step 2:** Connect the Home Medical Hub to the wireless network (Wi-Fi).

**WARNING:** Place the Home Medical Hub in a safe location so as not to fall on anyone and where the alarm sounds can be heard.

1. Following the on-screen instructions, connect the Home Medical Hub to a power outlet using the AC power cord and adapter and select **Continue**.

**Note:** The Home Medical Hub Status LED is white when powered on.

2. Hold your smart phone near the Home Medical Hub during setup. When complete, the App screen will display: *Success*.

**Note:** If the Bluetooth connection fails, hold the smart phone closer to the Home Medical Hub and try to reconnect. If the Home Medical Hub cannot connect, see ***Troubleshooting*** on page 39.

3. Next, connect to your Wi-Fi by selecting **Set Up Wi-Fi** and follow the app instructions.
4. Select the wireless *Network* and enter the *Password*. Select **Continue**.

**Note:** To view available Wi-Fi networks within range of the Home Medical Hub, click the *Network* field and select from the displayed list

5. Once the Wi-Fi connection is established, the *Home Medical Hub Wi-Fi Setup Successful* screen will appear. If the Home Medical Hub cannot connect, follow the on-screen instructions to verify the network and try again or view ***Troubleshooting*** on page 39.
6. Next, you will need to attach the Masimo Sensor and connect it with the Home Medical Hub. See ***Masimo Sensor Setup*** on page 22 for additional instructions.

## Masimo Sensor Setup

After the Home Medical Hub is set up, follow the on-screen instructions to:

- Set up the Masimo Sensor.
- Attach the sensor to your hand.
- Pair the sensor to the Home Medical Hub.
- Insert the Masimo Chip into the sensor.

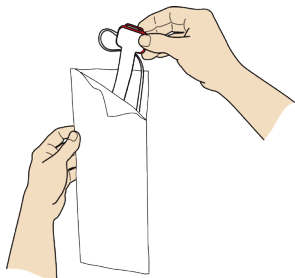
## Attaching the Masimo Sensor

Follow the instructions to attach the Masimo Sensor to your hand.


1. Gather all components required for monitoring. See ***Masimo Sensor and Masimo Chip Overview*** on page 18.


2. Remove the sensor from the pouch.

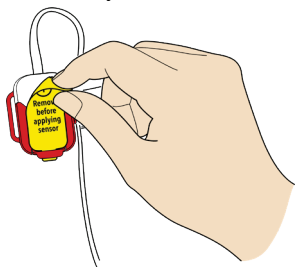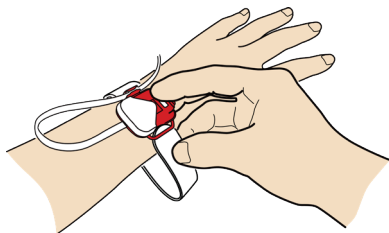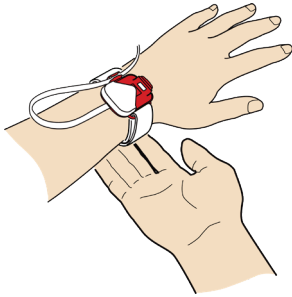
3. Peel off the yellow label to activate the battery.

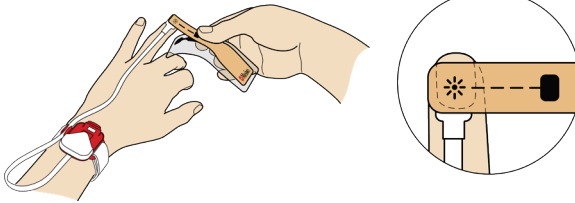
4. Place the sensor on the wrist as shown and carefully thread the tip of the strap through the open red c loop.


5. Wrap the loose attachment strap around the wrist and press to secure.


6. Remove part of the plastic film and place the star  on top of the fingernail of your ring finger.

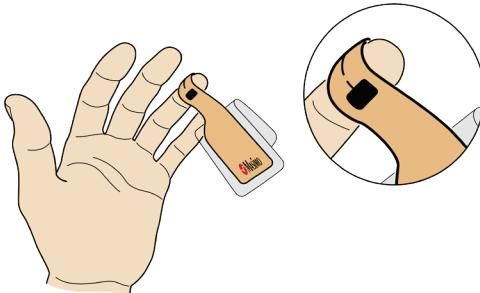
7. Wrap the tape around your finger so the square marking on the sensor tape is on the finger pad directly across from the star.

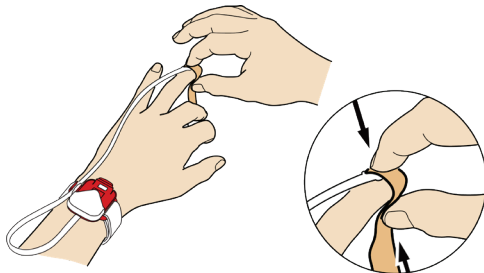
8. Squeeze top and bottom of finger to secure tape.

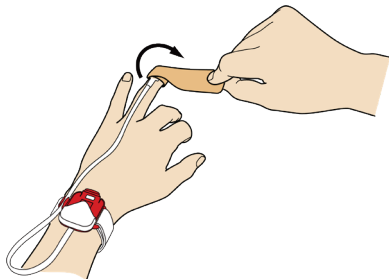9. Remove plastic film and continue to wrap remaining tape around your finger.

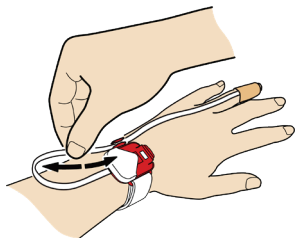
10. Adjust the sensor cable to a comfortable length. Go to ***Pairing the Masimo Sensor*** on page 26 for additional instructions.


## Pairing the Masimo Sensor

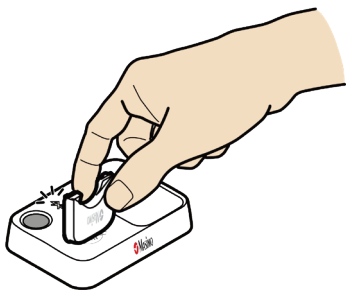
1. To pair the Masimo Sensor to the Home Medical Hub, hold the chip near the Home Medical Hub pairing symbol until an audible beep is heard and the Home Medical Hub Bluetooth LED flashes. See **Home Medical Hub Overview** on page 16.

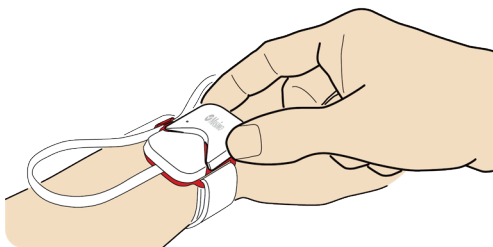2. Insert the Masimo Chip into open compartment on the sensor.

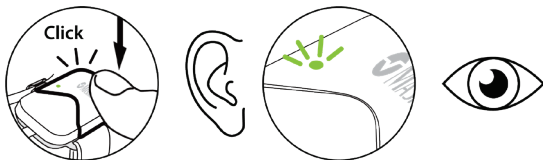3. Press down on the Masimo Chip to secure. Listen for a click and look for a green light on the chip to confirm successful chip connection to sensor.


4. Confirm successful pairing once the chip is inserted into the sensor by observing the light indicators.

#### **Home Medical Hub:**

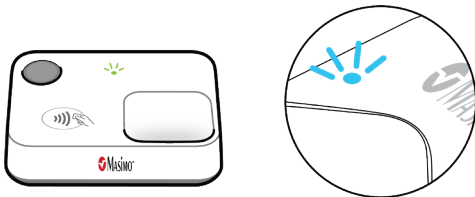- Bluetooth Status LED - turns off
- Status LED - turns green

#### **Masimo Chip:**

- Light Indicator - turns blue

- If the sensor and Home Medical Hub cannot pair, see **Troubleshooting** on page 39.

## Add Emergency Contacts

Follow the on-screen prompts to select your Emergency Contacts from your smart phone. This allows the Masimo SafetyNet Alert App to contact them in case of an emergency. An emergency contact is someone who will be alerted if you experience an Alert event.

1. Choose *Add From Contacts* or *Add Manually* to get started.
2. If prompted to *Allow Masimo SafetyNet Alert to access your contacts*, select **Allow**.  
**Note:** If you **Deny** access, Masimo SafetyNet Alert cannot add contacts for use in case of an emergency.
3. Select your emergency contacts from your contact list or add manually by entering the contact's name and mobile phone number.  
**CAUTION:** Do not use contacts with a landline phone number since they will not be able to receive text notifications in case of an emergency.
4. Select the types of alerts for the emergency contacts to receive and add the emergency contact.

Once contacts are added to receive alerts, the *Request Sent* screen appears. Select **Finish** to complete.

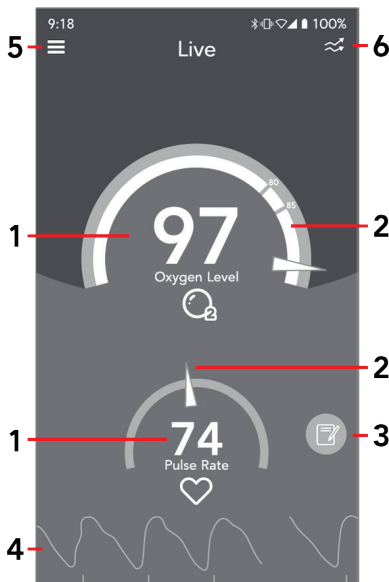
When finished, the *Live Data* screen appears. See **Step 3: View Live Data** on page 28.

- The emergency contact should receive a text message with a code. To accept, they must use the code and respond to the message as instructed.
- A text message is also sent to you when the emergency contact has accepted.

## Step 3: View Live Data

After setting up the Masimo SafetyNet Alert system, the app displays data from the Masimo Sensor. The *Live* screen is also the Masimo SafetyNet Alert App Main Screen, with access to other App functions.

\* If a sensor is not paired to the Home Medical Hub, a dash displays for the number, the pointer does not show on the gauge and the waveform does not display.

### 1. Data Display \*

Displays the data from the sensor. Touch the display to show additional information. See **Health Settings** on page 35.

### 2. Data Display Gauge \*

Displays the data from the sensor in a gauge style.

### 3. Timeline Notes Button

Touch to view and add Timeline notes.

### 4. Pulse Rate Waveform Display \*

Displays the waveform of the pulse rate. Touch the display to show additional information.

### 5. App Menu

Displays the *App Menu* and provides links to App functions and features.

### 6. Gauge/Trend View Toggle

Touch to toggle the screen to display data in gauge view (shown) or trend view.

# Alarms and Notifications

Masimo SafetyNet Alert provides visual and audible alarms when your oxygen level becomes low. Masimo SafetyNet Alert will also contact your emergency contacts and notify if your oxygen level drops too low and you are unaware or unable to respond.

## Alarm Levels and Messages

To help ensure your safety, a notification escalation policy is provided to establish three (3) levels of notifications. The level of notification is based on your risk level taking into account your oxygen level and length of time at a low oxygen level.

## Alert Events

Some Alert Events are pre-configured based upon your Option selection. For *COVID-19* and *Other* selections, Alert Event can be customized by you or on the advice of your physician. See ***Oxygen Level*** on page 36 for Alert Event defaults and available settings.

### Level 1

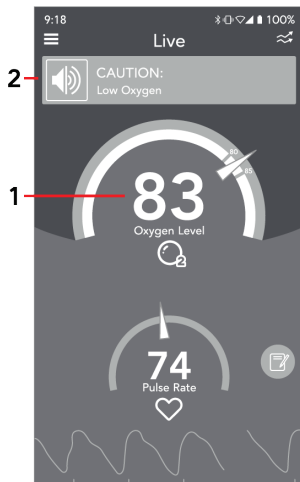
Visual Alarm	Audible Alarm	Emergency Contact
CAUTION (Yellow)	Yes	Not Contacted

When Level 1 is triggered, an alarm sounds on your Home Medical Hub and the App.

The App displays a yellow CAUTION banner at the top of screen.

Touch the Oxygen Level **(1)** or the yellow banner **(2)** to learn what the Level 1 notification means.

### LEVEL 1 CAUTION Alert Event Information

Information about the low oxygen level event is shown along with steps on how to help your oxygen level recovery.

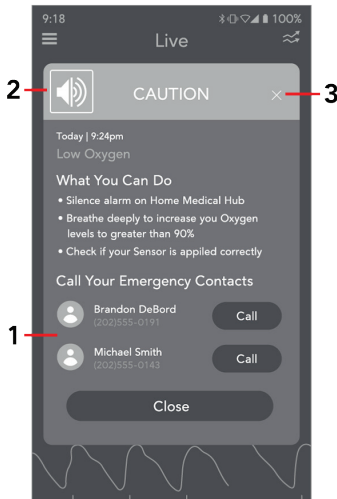
You have the option to:

1. Call a designated emergency contact **(1)** to assist you.
2. Silence the alarm

**On the App:** Tap the alarm icon on the top left corner of the pop-up **(2)**

**Home Medical Hub:** Press the button once.

See **Silence Alarms** on page 33 for additional information. To exit the screen, tap the "X" **(3)** or select the **Close** button.

### Levels 2 and 3

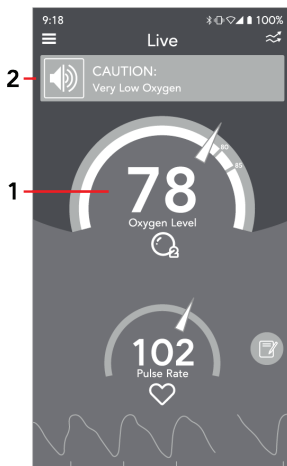
Notification Level	Visual Alarm	Audible Alarm	Emergency Contact
Level 2	WARNING (Red)	Yes	Message sent to Emergency Contact
Level 3	EMERGENCY (Flashing Red)		

When Level 2 or 3 alarms are triggered, an alarm will sound on your Home Medical Hub and the App.

The App will display your Oxygen Level in the red zone and a red banner at the top of screen.

- Level 2 Alarm: WARNING message.
- Level 3 Alarm: EMERGENCY message.

Touch the Oxygen Level **(1)** or red banner **(2)** to learn what the Level 2 or 3 alarm means.

## LEVEL 2 WARNING/LEVEL 3 EMERGENCY Alert Event Information

Information about the low oxygen level event is given along with steps on how to help your oxygen level recovery.

Designated emergency contacts are notified to assist you.

You have the option to:

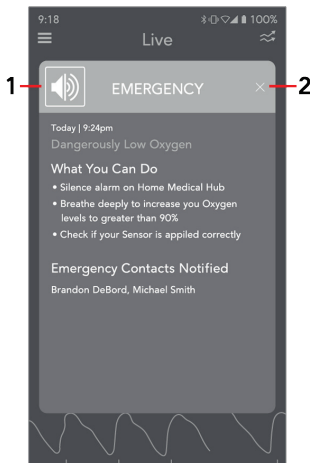
1. Silence the alarm

### Silence the Alarm

**On the App:** Tap the alarm icon on the top left corner of the pop-up (1)

**Home Medical Hub:** Press the button once.

See **Silence Alarms** on page 33 for additional information. To exit the screen, select the "X" (2) or select the **Close** button.

## Silence Alarms

An audio alarm can only be silenced for 2 minutes. After the 2 minutes, the alarm will sound again. You can continue to temporarily silence the alarm for 2 minutes at a time if needed. If your condition returns to a non-alarm level, the alarm will stop.


### Silencing Alarms From the App

Touch the Speaker icon  shown on the banner to silence the App alarm.

## Silencing Alarms From the Home Medical Hub

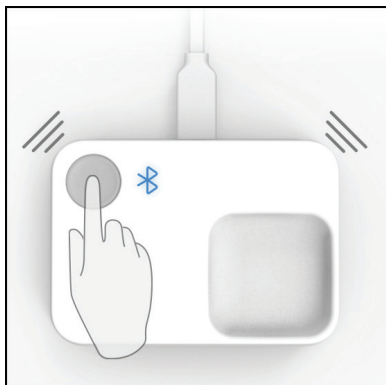Press and release the Home Medical Hub button to silence the Home Medical Hub alarm.

# Advanced Use

---

## Health Settings

The *Health Settings* menu is located within the *App Menu*. See **Step 3: View Live Data** on page 28. Options available include:

- Change Option - Information about available Options and how to change to different Option settings. See **Change Option** on page 35.
- Oxygen Level - Information about the Oxygen Level and available settings by Option. See **Oxygen Level** on page 36.
- Pulse Rate - Displays information about the Pulse Rate.
- Perfusion Index - Displays information about the Perfusion Index. The Perfusion Index on the Live Data screen can be turned On or Off from this screen.

## Change Option

When setting up the Masimo SafetyNet Alert App, you were asked to select an Option that best applied to you. The selected Option helps to identify the most relevant Alert Events for you at the time you selected it.

If your Option changes, you can perform the following:

1. Open the *App Menu*. See **Step 3: View Live Data** on page 28.
2. Select *Health Settings* from the *App Menu*.
3. In *Health Settings* select *Change Option*.
4. Under *Change Option*, select from the following conditions:
  - **COVID-19** - if you are recovering from COVID-19.
  - **Opioids** - if you are taking opioids.
  - **Other** - if you are not recovering from COVID-19 or not taking opioids.
5. Select OK.

The App will now update Alert Event options based upon your selection.

# Oxygen Level

The Oxygen Level screen shows the Oxygen Level (SpO<sub>2</sub>) Low Limit and the time below the Low Limit (Time Delay) for triggering the different Alert Events. The triggers are different depending on the Option selected.

Selecting the "?" icon next to the Option displays additional information about the Option.

## Low Oxygen Level and Alert Time Delay Default and Available Settings

For *COVID-19* or *Other*, you have the option to change the default Oxygen Level Low Limits and Alert Time Delays. Changes to the limit affect when Alert Events are triggered. Default settings can also be restored if needed.

### Alert Event Defaults and Settings (COVID-19)

Alert Type	Low SpO <sub>2</sub> Limit % Default	Low SpO <sub>2</sub> Limit % User Settings	Time Delay Default	Time Delay User Settings
Level 1 - Caution	93	Max: 95 Min: Warning Limit +1%	5 minutes	30 seconds to 5 minutes in 30 second increments
Level 2 - Warning	90	Max: Caution Limit –1% Min: Emergency Limit +1%	5 minutes	30 seconds to 5 minutes in 30 second increments
Level 3 - Emergency	80	Max: Warning Limit – 1% Min: 80	2 minutes	30 seconds to 2 minutes in 30 second increments

### Alert Event Defaults (Opioids)

Alert Type	Low SpO <sub>2</sub> Limit % Default	Low SpO <sub>2</sub> Limit % User Settings	Time Delay Default	Time Delay User Settings
Level 1 - Caution	85	NA	30 seconds	NA
Level 2 - Warning	80	NA	30 seconds	NA
Level 3 - Emergency	80	NA	180 seconds	NA
	60	NA	60 seconds	NA

### Alert Event Defaults (Other)

Alert Type	Low SpO <sub>2</sub> Limit % Default	Low SpO <sub>2</sub> Limit % User Settings	Time Delay Default	Time Delay User Settings
<b>Level 1 - Caution</b>	93	Max: 95 Min: Warning Limit + 1%	5 minutes	30 seconds to 5 minutes in 30 second increments
<b>Level 2 - Warning</b>	90	Max: Caution Limit - 1% Min: Emergency Limit + 1%	2 minutes	30 seconds to 5 minutes in 30 second increments
<b>Level 3 - Emergency</b>	80	Max: Warning Limit - 1% Min: 80	30 seconds	30 seconds to 2 minutes in 30 second increments

# Troubleshooting

## Masimo SafetyNet Alert Messages

The following section lists possible messages, the potential cause, and next steps.

Displayed Messages	Potential Causes	Next Steps
<i>"Difficulty in obtaining a reading"</i>	Interference while monitoring.	Please ensure that: <ol style="list-style-type: none"><li>1. The sensor is not on a finger with jewelry, such as a ring.</li><li>2. You move away from ambient lighting and computer displays/TVs.</li><li>3. You are not wearing artificial nails or excessive nail polish on the monitored finger.</li></ol> If you still experience issues, please contact Masimo Customer Support. See Customer Support.
<i>"Replace your sensor"</i>	Sensor is not working.	Replace sensor.
<i>"Wireless sensor disconnected during an alert"</i>	Sensor became disconnected during audible alerts for including: <ul style="list-style-type: none"><li>• Physiology alert or technical alert</li><li>• Low SpO<sub>2</sub></li><li>• Sensor off</li><li>• Obstructed battery</li><li>• Low battery</li><li>• Depleted battery</li></ul>	<ul style="list-style-type: none"><li>• Press the alarm silence button on the Home Medical Hub. Follow the instructions and attempt to pair the sensor again. See <b><i>Pairing the Masimo Sensor</i></b> on page 26.</li><li>• If you still experience issues, please contact Masimo Customer Support. See Customer Support.</li></ul>
<i>"Allow Masimo SafetyNet Alert to access your contacts"</i>	Access to contacts on smart phone has not been granted.	Allow smart phone permissions to share contacts with the Masimo SafetyNet Alert App.

Displayed Messages	Potential Causes	Next Steps
<i>"Your Emergency Contact list is empty. Touch the plus icon to start adding Emergency Contacts"</i>	When emergency contact list does not include a contact.	Follow instructions on the app to add emergency contacts.
<i>"Cannot modify a device that is in an active monitoring state. Please end your monitoring session and try again."</i>	When trying to change a setting that may disrupt alarms/notifications while in a monitoring state.	End the monitoring session before attempting to modify the device.
<i>"Confirm you Email"</i>	When an unverified email address was used to log into the app.	Confirm the email as requested.
<i>"A confirmation email was sent to emailaddress@example.com. Follow the instructions for confirming your email address then return to the App."</i>	When a new email address is submitted during setup.	Follow the instructions on the confirmation email.
<i>"An Emergency Contact is a designated friend, family member or caregiver that will be alerted in the event you need immediate assistance regarding your Oxygen levels. SafetyNet gives you the option to Add an Emergency Contact now or later."</i>	When adding an emergency contact to describe an emergency contact.	<ul style="list-style-type: none"> <li>• For information only, no action is required.</li> <li>• For more information on emergency contacts go to <b><a href="http://www.masimo.com/support">www.masimo.com/support</a></b>.</li> </ul>
<i>"Add Manually" "Create a new emergency contact. Please fill out the following information below"</i>	When adding an emergency contact manually.	For information only, no action is required.

Displayed Messages	Potential Causes	Next Steps
<p><i>"Request Sent"</i>  <i>"A request has been sent to your Emergency Contact. Follow up with your Emergency Contact to make sure they accept."</i></p>	<p>Displayed after successfully sending an emergency contact invitation</p>	<p>For information only, no action is required.</p>
<p><i>"Email Sent!"</i>  <i>"Please check your email to confirm"</i></p>	<p>Resend button was pressed to request a new email verification link.</p>	<p>For information only, no action is required.</p>
<p><i>"Device in use"</i>  <i>"Cannot modify a device that is actively monitoring. Please end your monitoring session and try again."</i></p>	<p>When an attempt was made to modify or delete the Home Medical Hub during active monitoring.</p>	<p>End the monitoring session before attempting to modify the device.</p>
<p><i>"Remove this emergency contact?"</i>  <i>"Remove selected emergency contacts?"</i></p>	<p>Displayed when removing one or more emergency contacts.</p>	<p>Select Remove to confirm or Cancel to cancel the operation.</p>
<p><i>"Remove this device?"</i></p>	<p>Displayed when attempting to remove Home Medical Hub</p>	<p>Select Remove to confirm or Cancel to cancel the operation.</p>
<p><i>"Save Changes?"</i>  <i>"Your unsaved changes will be lost. Save changes before closing?"</i></p>	<p>Displayed when attempting to update profile information.</p>	<p>Select Yes to confirm or No to not save changes.</p>
<p><i>"Your Home Medical Hub is paired"</i></p>	<p>The Home Medical Hub has paired successfully.</p>	<p>For information only, no action is required.</p>
<p><i>"No Home Medical Hub added"</i></p>	<p>Home Medical Hub is not connected to the app.</p>	<p>Follow instruction on the app or this manual to complete set up. See <b>Home Medical Hub Setup</b> on page 21.</p>

<b>Displayed Messages</b>	<b>Potential Causes</b>	<b>Next Steps</b>
<i>"Home Medical Hub not found"</i>	Timeout has occurred during setup when searching for Home Medical Hub.	Follow the troubleshooting section of this manual to resolve the issue.
<i>"Home Medical Hub Disconnected from Server"</i>	Home Medical Hub is disconnected from the cloud server.	Follow the troubleshooting section of this manual to resolve the issue.
<i>"Home Medical Hub Error"</i>	Home Medical Hub internal hardware fault.	The Home Medical Hub requires replacement, please contact Masimo Customer Support. See Customer Support.
<i>"Unable to connect to Wi-Fi"</i>	Wi-Fi connection failure due to timeout or invalid password.	Enter the correct password.
<i>"No Internet Connection"</i>	<ul style="list-style-type: none"> <li>• The smart phone is not connected to a Wi-Fi network.</li> <li>• The smart phone is not connected to a cellular network.</li> </ul>	Ensure the smart phone is connected to a Wi-Fi or cellular network.
<i>"Enable phone Bluetooth"</i>	Smart phone Bluetooth is turned off.	Turn on the smart phone Bluetooth.
<i>"The username entered already exists, please try another"</i>	User name already exist with another user.	Please choose another user name.
<i>"Battery low Warning"</i>	Sensor battery is low.	Replace sensor.
<i>"Depleted Battery"</i>	Sensor battery is depleted.	Replace sensor.
<i>"Wireless Sensor Disconnected"</i>	Sensor is not monitoring.	Ensure sensor placement steps are completed and you are obtaining readings.
<i>"Your oxygen level measurement has a low value of 76%"</i>	Level 1 alarm is triggered.	Follow instruction on the app to resolve the condition.

Displayed Messages	Potential Causes	Next Steps
<i>"Obstructed Battery Port"</i>	Sensor is not monitoring.	Please keep the sensor battery clear of any direct contact.
<i>"Place Sensor on Properly"</i>	Sensor is not monitoring.	Place the sensor on properly. You may go to the options menu to view the sensor placement videos. See <b>Masimo Sensor Setup</b> on page 22.
<i>"Allow Masimo SafetyNet Alert to access this device's location"</i>	Location services permissions have not yet been granted when attempting to search for a Bluetooth device.	Allow the smart phone permission to share its location with the Masimo SafetyNet Alert App.
<i>"Please Enable App Notifications"</i>	Notifications are not enabled on the smart phone for the Masimo SafetyNet Alert App.	Enable notifications on the smart phone for the Masimo SafetyNet Alert App.
<i>"Low Disc Space Warning!"</i>	Low diskspace during active monitoring and writing session trend.	Clear diskspace on the smart phone. A minimum of 100MB free disc space is required on the smart phone for the Masimo SafetyNet Alert App to operate properly.
<i>"You cannot access this item at the moment."</i>	Error when connecting to the server.	<ul style="list-style-type: none"> <li>• Try again.</li> <li>• If you still experience issues, please contact Masimo Customer Support. See Customer Support.</li> </ul>
<i>"There was an error processing your request."</i>	Error when connecting to the server.	<ul style="list-style-type: none"> <li>• Try again.</li> <li>• If you still experience issues, please contact Masimo Customer Support. See Customer Support.</li> </ul>
<i>"This request is not understood."</i>	Error when connecting to the server.	<ul style="list-style-type: none"> <li>• Try again.</li> <li>• If you still experience issues, please contact Masimo Customer Support. See Customer Support.</li> </ul>
<i>"Username or password invalid, please try again."</i>	Error when connecting to the server. Displayed when invalid credentials are used during the sign in process	<ul style="list-style-type: none"> <li>• Reenter the Username and Password.</li> <li>• If you still experience issues, please contact Masimo Customer Support. See Customer Support.</li> </ul>

<b>Displayed Messages</b>	<b>Potential Causes</b>	<b>Next Steps</b>
<i>"The username entered already exists, please try another."</i>	Error when connecting to the server. Displayed during registration when the username/email is already in use by another user.	Try a new Username and Password.
<i>"This account has already been verified."</i>	When attempting to resend an email Verification link to an already verified account.	Continue with Masimo SafetyNet Alert App setup, as email has already been verified.
<i>"There is no account registered under this username."</i>	Error when connecting to the server. Displayed when a user attempts to sign in with a username that does not exist.	Reenter the Username and Password. If you still experience issues, please contact Masimo Customer Support. See Customer Support.
<i>"User has already registered a different device"</i>	Error adding a new device. When adding a device fails due to different existing device	Follow steps to remove current device in order to add new device.
<i>"This Home Medical Hub has already been registered to an account."</i>	Error adding a device. When attempting to register a device that has already been registered to another user	Remove the device from the original account. If you still experience issues, please contact Masimo Customer Support. See Customer Support.
<i>"Device deleted successfully"</i>	When attempting to remove a device that no longer exists for this user on the cloud.	For information only, no action is required.
<i>"Invitation already sent to this Emergency Contact."</i>	When attempting to add an emergency contact that has a pending/active invitation with the current user.	For information only, no action is required.
<i>"Emergency Contact deleted successfully"</i>	When attempting to remove an emergency contact that no longer exists for this user on the cloud.	For information only, no action is required.
<i>"Maximum number of</i>	When attempting to add an	For information only, no action is required.

Displayed Messages	Potential Causes	Next Steps
<i>emergency contacts reached."</i>	emergency contact when the limit has reached. Limit = 10	
<i>"Something unexpected happened. Please try again later."</i>	System error.	<ul style="list-style-type: none"> <li>• Try again.</li> <li>• If you still experience issues, please contact Masimo Customer Support. See Customer Support.</li> </ul>

## Troubleshooting Masimo SafetyNet Alert

The following section lists possible symptoms, the potential cause, and next steps.

Symptom	Potential Causes	Next Steps
<i>Masimo SafetyNet Alert App does not turn on</i>	<ul style="list-style-type: none"> <li>• Depleted smart phone battery.</li> <li>• Masimo SafetyNet Alert App needs to be updated.</li> <li>• Incompatible smart phone.</li> </ul>	<ul style="list-style-type: none"> <li>• Connect smart phone to battery charger and charge battery. Refer to the smart phone's Operator's Manual or Directions For Use.</li> <li>• Check for Masimo SafetyNet Alert App updates.</li> <li>• Check smart phone compatibility. See <b>Specifications</b> on page 51.</li> <li>• Contact Masimo Customer Support. See Customer Support.</li> </ul>
<i>Masimo SafetyNet Alert App turns off</i>	<ul style="list-style-type: none"> <li>• Depleted smart phone battery.</li> <li>• Masimo SafetyNet Alert App needs to be updated.</li> <li>• Incompatible smart phone.</li> </ul>	<ul style="list-style-type: none"> <li>• Connect smart phone to battery charger and charge battery. Refer to the smart phone's Operator's Manual or Directions For Use.</li> <li>• Check smart phone compatibility. See <b>Specifications</b> on page 51.</li> <li>• Check for Masimo SafetyNet Alert App updates.</li> <li>• Contact Masimo Customer Support. See Customer Support.</li> </ul>

Symptom	Potential Causes	Next Steps
<p><i>Masimo SafetyNet Alert App does not communicate with Home Medical Hub</i></p>	<ul style="list-style-type: none"> <li>• Home Medical Hub is not powered on.</li> <li>• Smart phone is not in close proximity to Home Medical Hub.</li> <li>• Bluetooth on the smart phone is not turned on and/or not correctly configured.</li> <li>• Smart phone does not support Bluetooth Low Energy (BLE).</li> </ul>	<ul style="list-style-type: none"> <li>• Check that Home Medical Hub is plugged into the AC power supply.</li> <li>• Ensure the smart phone is in close proximity with Home Medical Hub.</li> <li>• Ensure Bluetooth on the smart phone is turned on.</li> <li>• Check smart phone compatibility. See <b>Specifications</b> on page 51.</li> <li>• Update the smart phone software. Refer to the smart phone's Operator's Manual or Directions For Use.</li> <li>• Contact Masimo Customer Support. See Customer Support.</li> </ul>
<p><i>Home Medical Hub does not connect to Wi-Fi or Masimo Cloud</i></p>	<ul style="list-style-type: none"> <li>• Home Medical Hub is not plugged in.</li> <li>• Incorrect Wi-Fi network selected.</li> <li>• Incorrect Wi-Fi password in entered.</li> <li>• Wi-Fi network is not correctly configured.</li> <li>• Masimo Cloud may be down.</li> </ul>	<ul style="list-style-type: none"> <li>• Ensure Home Medical Hub is plugged in.</li> <li>• Ensure smart phone is connected to correct Wi-Fi network. See <b>Home Medical Hub Setup</b> on page 21.</li> <li>• Ensure correct Wi-Fi network is selected. See <b>Home Medical Hub Setup</b> on page 21.</li> <li>• Ensure correct Wi-Fi password is entered. See <b>Home Medical Hub Setup</b> on page 21.</li> <li>• Check that the wireless features are correctly configured. Refer to the smart phone's Operator's Manual or Directions For Use.</li> <li>• Check network settings and availability. You may need to call network provider for further assistance.</li> <li>• Contact Masimo Customer Support. See Customer Support.</li> </ul>

Symptom	Potential Causes	Next Steps
<p><i>Masimo Sensor does not pair with Home Medical Hub</i></p>	<ul style="list-style-type: none"> <li>• Home Medical Hub is not plugged in.</li> <li>• Chip is not inserted into sensor.</li> <li>• Sensor is not in close proximity with Home Medical Hub during pairing.</li> <li>• Pairing button is released too soon.</li> <li>• Incorrect user logged into the Masimo SafetyNet Alert App.</li> <li>• Home Medical Hub has been registered with another account.</li> <li>• Depleted sensor battery.</li> </ul>	<ul style="list-style-type: none"> <li>• Ensure Home Medical Hub is plugged in.</li> <li>• Ensure the chip is firmly inserted into the sensor. See <b><i>Attaching the Masimo Sensor</i></b> on page 22.</li> <li>• Ensure the sensor is in close proximity with Home Medical Hub during pairing.</li> <li>• Ensure to press and hold the pairing button for 5 seconds until an audible tone is heard. See <b><i>Home Medical Hub Setup</i></b> on page 21.</li> <li>• Ensure correct user is logged into the Masimo SafetyNet Alert App.</li> <li>• Contact Masimo Customer Support. See Customer Support.</li> </ul>
<p><i>Home Medical Hub turns off (not lit up)</i></p>	<ul style="list-style-type: none"> <li>• Home Medical Hub is not plugged in.</li> <li>• Internal components may not be working properly.</li> </ul>	<ul style="list-style-type: none"> <li>• Ensure Home Medical Hub is plugged in.</li> <li>• Contact Masimo Customer Support.</li> </ul>
<p><i>Home Medical Hub speaker does not work (no sound or muffled sound)</i></p>	<ul style="list-style-type: none"> <li>• Speaker may be blocked by environment (ex. blankets or other appliances).</li> <li>• Internal components may not be working properly.</li> </ul>	<ul style="list-style-type: none"> <li>• Turn Home Medical Hub on and off by unplugging the unit.</li> <li>• Check that the device speaker is not being muffled.</li> <li>• Check that the Home Medical Hub is on a flat surface with minimal objects surrounding it.</li> <li>• Contact Masimo Customer Support. See Customer Support.</li> </ul>

<b>Symptom</b>	<b>Potential Causes</b>	<b>Next Steps</b>
<i>Masimo Cloud unavailable</i>	<ul style="list-style-type: none"> <li>• Wi-Fi is not turned on and/or not correctly configured.</li> <li>• Wireless service is weak or unavailable in the current location.</li> </ul>	<ul style="list-style-type: none"> <li>• Ensure the smart phone is within range of the wireless network for connection to the system.</li> <li>• Check that the wireless feature for smart phone is on and correctly configured. Refer to the smart phone's Operator's Manual or Directions For Use.</li> <li>• Check Wi-Fi network settings and availability.</li> <li>• Check wireless availability for location.</li> <li>• Update the smart phone software. Refer to the smart phone's Operator's Manual or Directions For Use.</li> <li>• Contact Masimo Customer Support. See Customer Support.</li> </ul>
<i>Incorrect user data is displayed</i>	<ul style="list-style-type: none"> <li>• Incorrect user currently logged into the Masimo SafetyNet Alert App.</li> </ul>	<ul style="list-style-type: none"> <li>• Ensure the correct user is logged into the Masimo SafetyNet Alert App.</li> <li>• Log out the incorrect user and log into the correct account.</li> <li>• Restart Masimo SafetyNet Alert App and login to the system.</li> <li>• Contact Masimo Customer Support. See Customer Support.</li> </ul>
<i>NO user data is displayed</i>	<ul style="list-style-type: none"> <li>• Incorrect user logged into the Masimo SafetyNet Alert App.</li> <li>• Wi-Fi is not correctly configured.</li> <li>• Smart phone settings are incorrect.</li> <li>• Sensor is not connected to the Home Medical Hub.</li> <li>• No previous sessions have been recorded.</li> </ul>	<ul style="list-style-type: none"> <li>• Ensure the correct user is logged into the Masimo SafetyNet Alert App.</li> <li>• Restart Masimo SafetyNet Alert App and login to the system.</li> <li>• Check that the wireless feature is correctly configured. Refer to the smart phone's Operator's Manual or Directions For Use.</li> <li>• Check network settings and availability. May have to call network provider for further assistance.</li> <li>• Contact Masimo Customer Support. See Customer Support.</li> </ul>

Symptom	Potential Causes	Next Steps
<i>Delayed data updates</i>	<ul style="list-style-type: none"> <li>• Wi-Fi is not correctly configured.</li> </ul>	<ul style="list-style-type: none"> <li>• Check that the wireless feature is correctly configured. Refer to the smart phone's Operator's Manual or Directions For Use.</li> <li>• Check network settings and availability. You may need to call network provider for further assistance.</li> <li>• Restart Masimo SafetyNet Alert App and login to the system. See</li> <li>• Contact Masimo Customer Support. See Customer Support.</li> </ul>
<i>User alarms do not appear during events</i>	<ul style="list-style-type: none"> <li>• Incorrect user currently logged into Masimo SafetyNet Alert App.</li> <li>• Wi-Fi network not available.</li> <li>• System settings have changed.</li> </ul>	<ul style="list-style-type: none"> <li>• Ensure Masimo SafetyNet Alert is connected to Wi-Fi. See <b>Home Medical Hub Setup</b> on page 21.</li> <li>• Ensure the correct user is logged into Masimo SafetyNet Alert App.</li> <li>• Restart Masimo SafetyNet Alert App and login to the system.</li> <li>• Ensure alert notification settings are turned on.</li> <li>• Contact Masimo Customer Support. See Customer Support.</li> </ul>
<i>Masimo SafetyNet Alert App does not detect that the sensor is applied to the hand</i>	<ul style="list-style-type: none"> <li>• Sensor not properly placed on user.</li> <li>• Sensor not properly paired to Home Medical Hub.</li> <li>• Damaged sensor.</li> <li>• Internal failure.</li> </ul>	<ul style="list-style-type: none"> <li>• Reapply the sensor.</li> <li>• Pair the sensor with Home Medical Hub. See <b>Pairing the Masimo Sensor</b> on page 26.</li> <li>• Replace sensor.</li> <li>• Turn Home Medical Hub Off and On by unplugging the device then plug it in to start up</li> <li>• Contact Masimo Customer Support. See Customer Support.</li> </ul>

## Troubleshooting Measurements

The following section lists possible measurement symptoms, potential causes, and next steps. For more information, see **Safety Information** on page 7.

Symptom	Potential Causes	Next Steps
<i>Difficulty obtaining a reading.</i>	<ul style="list-style-type: none"><li>• Incorrect placement of sensor on user.</li><li>• Misalignment of sensor components.</li><li>• Low perfusion (blood flow).</li><li>• Excessive user motion.</li><li>• Excessive ambient or strobing light.</li><li>• Low battery/SafetyNet Home Medical Hub not plugged into AC power supply.</li></ul>	<ul style="list-style-type: none"><li>• Check the placement and alignment of the sensor on the hand. Re-apply sensor or move to a different location.</li><li>• Allow time for the parameter measurement to stabilize.</li><li>• Check if blood flow to the sensor location is restricted. Warm the hand where the sensor is placed.</li><li>• Minimize or eliminate motion at the monitoring location.</li><li>• Shield the sensor from excessive or strobing light.</li><li>• Replace sensor.</li><li>• Contact Masimo Customer Support. See Customer Support.</li></ul>
<i>Measurement values displayed as dashes.</i>	<ul style="list-style-type: none"><li>• Measurement may still be in progress.</li><li>• Incorrect placement of sensor on user.</li><li>• No Connection to Masimo Cloud.</li><li>• Sensor is damaged, not functioning or has a dead battery.</li></ul>	<ul style="list-style-type: none"><li>• Allow time for the parameter measurement to stabilize.</li><li>• Check the placement of the sensor on the hand. Re-apply sensor or move to a different location.</li><li>• Check if blood flow to the sensor location is restricted.</li><li>• Replace sensor.</li><li>• Contact Masimo Customer Support. See Customer Support.</li></ul>
<i>Unexpected or unlikely measurement values</i>	<ul style="list-style-type: none"><li>• Incorrect placement of sensor on user.</li><li>• Low signal quality.</li></ul>	<ul style="list-style-type: none"><li>• Check the placement of the sensor on the hand. Re-apply sensor or move to a different location.</li><li>• Move the sensor to a location on the user's body with stronger blood flow. Warm the hand where the sensor is placed.</li><li>• Contact Masimo Customer Support. See Customer Support.</li></ul>

# Appendix

---

## Specifications

### Masimo SafetyNet Alert App

#### Measurement Range

Measurement	Display Range	Unit of Measure
Oxygen Level (SpO <sub>2</sub> )	0 to 100	%
Pulse Rate (PR)	0 to 240	bpm
Perfusion Index (PI)	0 to 20	N/A

#### Smart Phone Compatibility

Item	Specification*
Operating System	Android 6.0 (Marshmallow) (minimum)
	iOS 12.0 (minimum)

\* For complete specifications, refer to [www.masimo.ca/support](http://www.masimo.ca/support).

## Masimo Sensor

### Weight Range and Measurement Site

Population	Weight	Measurement Site
Adult and Pediatric	> 40 kg (> 88.2 lbs)	Finger
Neonate	< 3 kg (< 6.6lbs)	Finger or Foot

## Accuracy (ARMS\*)

Oxygen Level (SpO <sub>2</sub> )		
Range	70% to 100%	
No Motion [1]	Adults, Pediatrics	2%
	Neonates	3%
Motion [2]	Adults, Pediatrics, Neonates	3%
Low perfusion [3]	Adults, Pediatrics	2%
	Neonates	3%
Pulse Rate (PR)		
Range	25 bpm to 240 bpm	
No motion	Adults, Pediatrics, Neonates	3 bpm
Motion	Adults, Pediatrics, Neonates	5 bpm
Low Perfusion [4]	Adults, Pediatrics, Neonates	3 bpm

\*  $A_{RMS}$  accuracy is a statistical calculation of the difference between device measurements and reference measurements. Approximately two-thirds of the device measurements fell within  $\pm A_{RMS}$  of the reference measurements in a controlled study.

## Oxygen Level (SpO<sub>2</sub>) ARMS Performance Specifications

The tables below provides  $A_{RMS}$  (Accuracy Root Mean Square) values measured using the Masimo Sensor (Radius PPG) under no motion, with Masimo Technology in a clinical study.

Measurement $A_{RMS}$ Values for Radius PPG Sensors	
SpO <sub>2</sub> Accuracy Range (%)	ARMS (%)
90-100	1.14
80-90	1.29
70-80	1.41

Measurement $A_{RMS}$ Values for Radius PPG Sensors	
SpO <sub>2</sub> Accuracy Range (%)	ARMS (%)
70-100	1.33

The below Bland-Altman plot represents the correlation of the  $(SpO_2 + SaO_2)/2$  versus  $(SpO_2 - SaO_2)$  under no motion with an upper 95% and lower 95% limits of agreement.

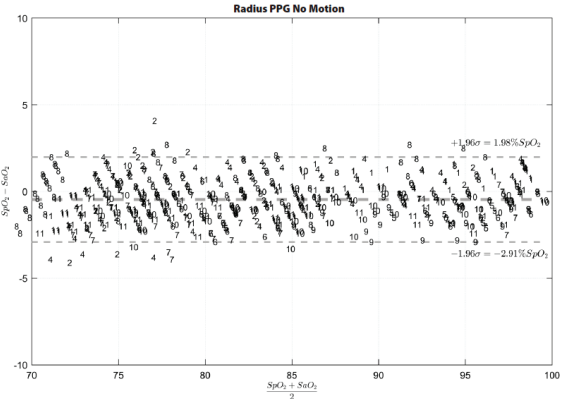
Figure 1: Radius PPG Sensors (ARMS 70-100%)

### Electrical

Battery - Masimo Sensor	
Run Time	96 hours in typical continuous usage

### Environmental

Masimo Sensor Environmental Conditions	
Operating Temperature	32°F to 104°F (0°C to 40°C)

<b>Masimo Sensor Environmental Conditions</b>	
Storage Temperature	32°F to 122°F (0°C to 50°C)
Operating Humidity	5% to 95%, non-condensing
Storage Humidity	5% to 95%, non-condensing

## Masimo Chip

### Environmental

<b>Masimo Chip Environmental Conditions</b>	
Operating Temperature	32°F to 104°F (0°C to 40°C)
Storage Temperature	-40°F to 158°F (-40°C to 70°C)
Operating Humidity	5% to 95%, non-condensing
Storage Humidity	5% to 95%, non-condensing
Atmospheric Pressure	540 to 1060 mBar @ ambient temperature and humidity

## Home Medical Hub

### Electrical

<b>AC Power Requirements</b>	
AC Power Input (External Power Supply)	100 to 240 VAC, 50 to 60 Hz, 1.2 A
DC Power Input (Home Medical Hub)	5 VDC, 750 mA

## Physical Characteristics

<b>Home Medical Hub</b>	
Dimensions	3.45" x 2.48" x 0.78" (8.76 cm x 6.3 cm x 2 cm)

<b>Home Medical Hub</b>	
Weight	0.12 lbs. (56g)

## Environmental

<b>Environmental Conditions - Home Medical Hub</b>	
Operating Temperature	41°F to 104°F (5°C to 40°C)
Storage Temperature	-13°F to 158°F (-25°C to 70°C)
Operating Humidity	10% to 95%, non-condensing
Storage Humidity	10% to 95%, non-condensing

## Compliance

<b>EMC Compliance</b>
IEC 60601-1-2:2014, Class B

<b>Safety Standards Compliance</b>
IEC 60601-1
IEC 62304
IEC 60601-1-11

<b>Equipment Classification per IEC 60601-1</b>	
Type of Protection	Class II (AC Power)
Degree of Protection against Electrical Shock	Type BF-Applied Part
Protection against harm from Water and Particulate Matter	IP22 (Protection from solid foreign objects $\geq 12.5$ mm diameter and against ingress from vertically falling water drops when enclosure tilted up to 15°)

**Equipment Classification per IEC 60601-1**

Mode of Operation

Continuous

## Wireless Specifications

**Communication (Bluetooth)**

Type	Bluetooth
Frequency	2402-2480 MHz
Max Peak Output Power	8.45 dBm
Classification of Output Power Rating	Conducted
Output Power Type	Fixed at the Factory
Modulation Types	GFSK
Modulation Signals	Analog and Digital
Available Data Rates	1 Mbps
Recommended Max. Range	100 ft (~30 meters) line-of-sight

**Communication (Wi-Fi)**

Type	WLAN Radio: IEEE 802.11 b/g/n
Frequency	802.11b/g/n(HT20): 2412-2462 MHz 802.11n(HT40): 2422-2452 MHz
Max Peak Output Power	WLAN 27.12 dBm
Classification of Output Power Rating	Conducted
Output Power Type	Fixed at the Factory
Modulation Types	802.11b: DSSS 802.11g/n(HT20/HT40): OFDM
Modulation Signals	Analog and Digital

**Communication (Wi-Fi)**

Available Data Rates	802.11b - 1, 2, 5.5, 11 Mbps. 802.11g - 6, 9, 12, 18, 24, 36, 48, 54 Mbps 802.11n- MCS0 – MCS7
----------------------	--

**Security and Authentication**

Encryption	64/128-bit WEP, Dynamic WEP, WPA-TKIP, WPA2-AES
Authentication	Open System, Shared Key, Pre-Shared Key (PSK), 802.1X: LEAP, PEAP, TTLS, TLS, EAP-FAST

**Radio Compliance**

USA	Contains FCC ID: 2AC7Z-ESP32SOLO1 FCC ID: VKF-CONNHUB
Canada	Contains IC: 21098-ESP32SOLO1 IC: 7362A-CONNHUB
Europe	EU Radio Equipment Directive (RED 2014/53/EU) EN 300 330 V2.1.1:2017 EN 301 489-3 V2.1.1:2019 1999/519/EC EN 62311:2020

# Guidance and Manufacturer's Declarations - Electromagnetic Compliance

## Electromagnetic Emissions

<b>Electromagnetic Emissions</b>		
The ME Equipment is intended for use in the electromagnetic environment specified below. The customer or the user of the ME Equipment should assure that it is used in such an environment.		
<b>Emission Test</b>	<b>Compliance</b>	<b>Electromagnetic Environment - Guidance</b>
RF Emissions (Radiated) CISPR 11	Group 1 Class B	ME Equipment uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF Emissions (Conducted) CISPR 11	Group 1 Class B	Suitable for use in all establishments, including domestic environments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.
Harmonic Emissions IEC 61000-3-2	Class A	
Voltage fluctuations/ Flicker emissions IEC 61000-3-3	Complies	

## Electromagnetic Immunity

<b>Electromagnetic Immunity</b>			
The ME Equipment is intended for use in the electromagnetic environment specified below. The customer or the user of the ME Equipment should assure that it is used in such an environment.			
<b>Immunity Test</b>	<b>IEC 60601 Test Level</b>	<b>Compliance Level</b>	<b>Electromagnetic Environment - Guidance</b>
Electrostatic discharge (ESD) IEC 61000-4-2	+/- 8 kV contact +/- 15 kV air	+/- 8 kV contact +/- 15 air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%.
Electrical fast transient/ burst IEC 61000-4-4	+/- 2 kV for power lines +/- 1 kV for input/ output lines	+/- 2 kV for power lines +/- 1 kV for input/ output lines	Mains power quality should be that of a typical commercial or hospital environment.
Surge IEC 61000-4-5	+/- 1 kV line(s) to line(s) +/- 2 kV line(s) to earth	+/- 1 kV line(s) to line(s) +/- 2 kV line(s) to earth	Mains power quality should be that of a typical commercial or hospital environment.
Conducted RF IEC 61000-4-6	3 Vrms	3 Vrms	Performed over 0.15-80 MHz
	6 Vrms	6 Vrms	Performed on the following ISM (industrial, scientific and medical) bands of frequency: The bands between 0,15 MHz and 80 MHz are 6,765 MHz to 6,795 MHz; 13,553 MHz to 13,567 MHz; 26,957 MHz to 27,283 MHz; and 40,66 MHz to 40,70 MHz. The amateur radio bands between 0,15 MHz and 80 MHz are 1,8 MHz to 2,0 MHz, 3,5 MHz to 4,0 MHz, 5,3 MHz to 5,4 MHz, 7 MHz to 7,3 MHz, 10,1 MHz to 10,15 MHz, 14 MHz to 14,2 MHz, 18,07 MHz to 18,17 MHz, 21,0 MHz to 21,4 MHz, 24,89 MHz to 24,99 MHz, 28,0 MHz to 29,7 MHz and 50,0 MHz to 54,0 MHz

<b>Electromagnetic Immunity</b>			
Power frequency (50 / 60 Hz) magnetic field IEC 61000-4-8	30 A/m	30 A/m	Power frequency magnetic fields should be at levels characteristic of typical location in a typical hospital environment.
Voltage dips on power supply input lines IEC 61000-4-11	0% $U_T^1$ , 0.5 cycle, at 0°, 45°, 90°, 135°, 180°, 225°, 270°, and 315°; 0% $U_T$ 1 cycle, and 70% $U_T$ 25/30 cycles at 0°	0% $U_T^1$ , 0.5 cycle, at 0°, 45°, 90°, 135°, 180°, 225°, 270°, and 315°; 0% $U_T$ 1 cycle, and 70% $U_T$ 25/30 cycles at 0°	Mains power quality should be that of a typical commercial or hospital environment.
Voltage Interruptions on power supply input lines IEC 61000-4-11	0% $U_T$ , 250/300 cycle	0% $U_T$ , 250/300 cycle	
Radiated RF IEC 61000-4-3	10 V/m	10 V/m	Performed over 80 MHz to 2.7 GHz
<p><b>Note 1:</b> At 80 MHz and 800 MHz, the higher frequency range applies.</p> <p><b>Note 2:</b> These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.</p>			
<p>Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the ME Equipment is used exceeds the applicable RF compliance level above, the ME Equipment should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as re-orienting or relocating the ME Equipment.</p>			

<sup>1</sup>  $U_T$ : Rated voltage for the equipment.

## Test Specifications for ENCLOSURE PORT IMMUNITY to RF Wireless Communication Equipment

Test Frequency (MHz)	Band (a) (MHz)	Service (a)	Modulation (b)	Maximum Power (W)	Distance (m)	Immunity Test Level (V/m)
385	380-395	TETRA 400	Pulse modulation (b) 18 Hz	1,8	0,3	27
450	430-470	GMRS 460, FRS 460	FM (c) +/- 5 kHz deviation 1 kHz sine	2	0,3	28
710	704-787	LTE Band 13, 17	Pulse modulation (b) 217 Hz	0,2	0,3	9
745						
780						
810	800-960	GSM 800/900, TETRA 800, iDEN 820, CDMA 850, LTE Band 5	Pulse modulation (b) 18 Hz	2	0,3	28
870						
930						
1 720	1 700-1 990	GSM 1800; CDMA 1900; GSM 1900; DECT; LTE Band 1, 3. 4. 35: UMTS	Pulse modulation (b) 217 Hz	2	0,3	28
1 845						
1 970						
2 450	2 400-2 570	Bluetooth, WLAN, 802.11 b/g/n, RFID 2450, LTE Band 7	Pulse modulation (b) 217 Hz	2	0,3	28
5 240	5 100-5 800	WLAN 802.11 a/n	Pulse modulation (b) 217 Hz	0,2	0,3	9
5 500						
5 785						

Test Frequency (MHz)	Band (a) (MHz)	Service (a)	Modulation (b)	Maximum Power (W)	Distance (m)	Immunity Test Level (V/m)
----------------------	----------------	-------------	----------------	-------------------	--------------	---------------------------

**Note:** If necessary to achieve the IMMUNITY TEST LEVEL, the distance between the transmitting antenna and the ME EQUIPMENT or ME SYSTEM may be reduced to 1 m. The 1 m test distance is permitted by IEC 61000-4-3.

- (a) For some services, only the uplink frequencies are included.  
 (b) The carrier shall be modulated use a 50% duty cycle square wave signal.  
 (c) As an alternative to FM modulation, 50% pulse modulation at 18 Hz may be used because while it does not represent actual modulation, it would be worst case.

## Recommended Separation Distances

### Recommended Separation Distance Between Portable and Mobile RF Communication Equipment and the ME Equipment

The ME Equipment is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the ME Equipment can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the ME Equipment as recommended below, according to the maximum output power of the communication equipment.

Rated maximum output power of transmitter (W)	Separation Distance According to Frequency of Transmitter (m)
	$d = 0.6 * \text{Sqrt}(P)$
0.01	0.06
0.1	0.19
1	0.6
10	1.9
100	6

For transmitters rated at a maximum output power not listed above, the recommended separation distance  $d$  in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where  $P$  is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.

## Recommended Separation Distance Between Portable and Mobile RF Communication Equipment and the ME Equipment

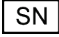
**Note 1:** At 80 MHz and 800 MHz, the higher frequency range applies.

**Note 2:** These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

## Symbols

The following symbols may appear on the product or product labeling:

Symbol	Description	Symbol	Description
	Follow instructions for use		Consult instructions for use
	Mark of conformity to European medical device directive 93/42/EEC		Separate collection for electrical and electronic equipment (WEEE)
IP22	Protection from solid foreign objects $\geq 12.5$ mm diameter and against ingress from vertically falling water drops when enclosure tilted up to $15^\circ$	IP47	Protection from solid foreign objects $< 1$ mm diameter and against immersion between 15 centimeters and 1 meter in depth
	Non-Sterile		Federal Communications Commission (FCC) Licensing
	Authorized representative in the European community		Recyclable
IC Model:	Innovation, Science and Economic Development Canada (ISED)	FCC ID:	Identifies unit has been registered as a radio device
	Caution		Product contains no PVC (polyvinyl chloride) material
	Manufacturer		Not made with natural rubber latex

Symbol	Description	Symbol	Description
	Date of manufacture YYYY-MM-DD		Catalog number (model number)
	Storage temperature range		Masimo reference number
	Keep dry		Serial number
	Storage humidity limitation		Do not use if package is damaged
	Atmospheric pressure limitation		AC current
	Wireless Symbol level		The names and content of the toxic and hazardous substances or elements shall be provided in the product instruction manual
	China Restriction of Hazardous Substances	--	--
	Instructions/Directions for Use/Manuals are available in electronic format @ <a href="http://www.Masimo.com/TechDocs">http://www.Masimo.com/TechDocs</a> <b>Note:</b> eFU is not available in all countries.		

## Citations

[1] The Masimo SET Technology has been validated for no motion accuracy in human blood studies on healthy adult male and female volunteers with light to dark pigmented skin in induced hypoxia studies in the range of 70%-100% SpO<sub>2</sub> against a laboratory co-oximeter.

[2] The Masimo SET Technology has been validated for motion accuracy in human blood studies on healthy adult male and female volunteers with light to dark pigmented skin in induced hypoxia studies while performing rubbing and tapping motions, at 2 to 4 Hz at an amplitude of 1 to 2 cm and a non-repetitive motion between 1 to 5 Hz at an amplitude of 2 to 3 cm in induced hypoxia studies in the range of 70%-100% SpO<sub>2</sub> against a laboratory co-oximeter.

[3] The Masimo SET Technology has been validated for low perfusion accuracy in bench top testing against a Biotek Index 2 simulator and Masimo's simulator with signal strengths of greater than 0.02% and transmission of greater than 5% for saturations ranging from 70% to 100%.

[4] The Masimo SET Technology has been validated for pulse rate accuracy for the range of 25-240 bpm in bench top testing against a Biotek Index 2 simulator and Masimo's simulator with signal strengths of greater than 0.02% and transmission of greater than 5% for saturations ranging from 70% to 100%.

\*Registered trademark of Fluke Biomedical Corporation, Everett, Washington.

## Service and Maintenance

### Cleaning

The Home Medical Hub and Masimo Chip are reusable devices. The devices are supplied and used non-sterile.

The Masimo Chip should be cleaned before and after it has been applied to a user and/or in accordance with local and governmental regulations to minimize the risk of cross-contamination.

### Smart Phone Cleaning

To properly clean the smart phone, refer to the smart phone's Operator's Manual or Directions For Use.

### Home Medical Hub Cleaning

**CAUTION:** Check the Home Medical Hub for possible cracks or opening before cleaning.

**CAUTION:** Do not allow liquids to enter the interior of the Home Medical Hub.

The outer surfaces can be cleaned either with a soft cloth dampened with a mild detergent and warm water solution or they can be wiped down with the following cleaning solutions:

- 70% isopropyl alcohol
- 1:10 bleach to water solution (0.5% sodium hypochlorite)
- Super Sani-Cloth® Wipes (55% isopropyl alcohol, 0.5% quaternary ammonium chloride)
- Windex® (1.5% 2-Butoxyethanol, 1.5% ethylene glycol hexyl ether, 5% isopropyl alcohol)

- Formula 409® Antibacterial All-Purpose Cleaner (1.5% Lauramine oxide, 0.4% n-Alkyl dimethyl benzyl ammonium chloride)

## Masimo Chip Cleaning

**WARNING:** Before cleaning, make sure the sensor and chip are not applied to the hand.

To surface clean the chip and wireless receiver:

1. Wipe all surfaces of the chip and wireless receiver with one of the following:
  - a. 70% Isopropyl alcohol
  - b. 10% (1:10) chlorine bleach to water solution
  - c. Quaternary ammonium chloride solution
2. Inspect for visible debris and repeat the above cleaning step as needed.
3. Dry cleaned parts before use.

### CAUTIONS:

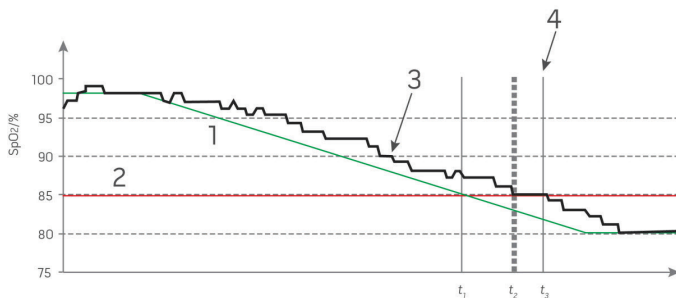
- To avoid permanent damage to the chip and wireless receiver, do not use undiluted bleach (5% - 5.25% sodium hypochlorite) or any other cleaning solution not recommended.
- Do not immerse the chip and wireless receiver in any liquid solution.
- Do not sterilize by irradiation, steam, autoclave or ethylene oxide.

## Customer Support

For answers to frequently asked questions (FAQ's) and product support, along with troubleshooting for your Masimo SafetyNet Alert product, please go to [www.masimo.ca/support](http://www.masimo.ca/support) or email [customer@masimo.com](mailto:customer@masimo.com).

## Concepts of Alarm Response Delay

As with any pulse oximeter equipment, the audible and visual alarms are subject to alarm response delay, which is composed of Alarm Condition Delay and Alarm Signal Generation Delay. Alarm Condition Delay is the time from the occurrence of the triggering event to when the alarm system determines the alarm condition exists. While Alarm Signal Generation Delay is the time from the onset of an alarm condition to the generation of its alarm signal. The graphic below is a simplified illustration of the concept of alarm response delay and does not reflect actual lengths of delays.

Reference	Definition	Reference	Definition
1	SaO <sub>2</sub>	4	Alarm Signal Generation
2	Alarm Limit	SpO <sub>2</sub>	Saturation
3	Displayed SpO <sub>2</sub>	t	Time

The Alarm Condition Delay is graphically represented as  $t_2 - t_1$  in the figure above to show the delay due to processing and averaging.

The Alarm Signal Generation Delay is graphically represented as  $t_3 - t_2$  in the figure above to show the delay due to alarm system strategy and communication time.

The overall alarm system delay time is graphically represented as  $t_3 - t_1$ .


For more information about alarm response delay, refer to ISO 80601-2-61.

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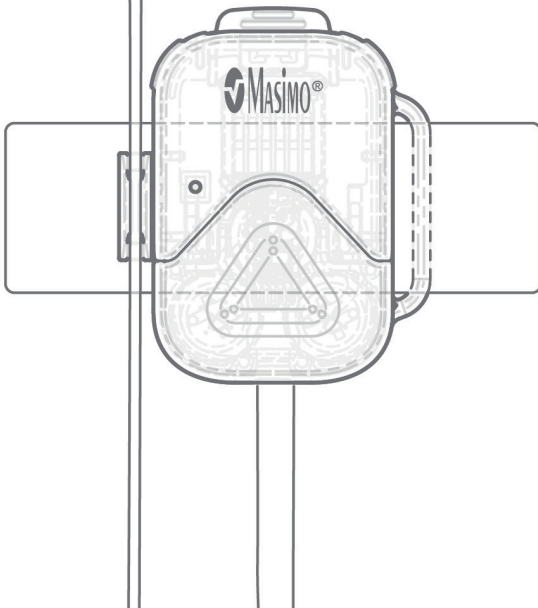## G

Getting Started • 19  
Getting Started Guides • 5  
Guidance and Manufacturer's Declarations  
- Electromagnetic Compliance • 58  
Guides de démarrage • 74

## H

Health Settings • 29, 35  
Home Medical Hub • 54  
Home Medical Hub Cleaning • 65  
Home Medical Hub Lights • 16, 17  
Home Medical Hub Overview • 16, 26  
Home Medical Hub Setup • 21, 41, 46, 47,  
49

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## S

- Safety Information • 7, 50
- Safety Warnings • 7, 76

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## W

Weight Range and Measurement Site • 51

Who the Device is intended to be used on  
(Indications for Use) • 6

Wireless Specifications • 56

## Y

Your Body and Oxygen • 12

[www.masimo.com](http://www.masimo.com)

303565/LAB-11582A-0322